

See Our Products in a
Whole New **Light!**



Spring 2011
Catalog

NEW! State-of-the-Art Mobile Dental Cart



AMC-20 Mobile Dental Cart

\$9,400

- Fully self-contained ergonomic dental cart with built-in air compressor
- 2 automatic highspeed air handpiece connections
- 1 automatic lowspeed handpiece connection
- High & low volume suction with automatic overflow protection and solids trap
- 3-way air/water syringe
- Unique pivoting control arms & easy grip side handles
- Two 1 liter bottles for irrigation solutions
- Large counter
- 2 duplex accessory hospital grade outlets
- Air filter moisture separator
- Integrated electric waste removal pump
- Handpiece flush system
- External quick connects for auxiliary air/water
- Resettable circuit breaker protection
- Sound insulated for super quiet operation
- Easy access for maintenance
- Wet/dry foot control
- Easy access to controls
- Large casters
- Weight: 120 lbs (54.4 Kg)



AMC-20 Options

• ASC-10-20CK	Piezo Scaler\$1,295
• AA-20CL-CK	Curing Light\$1,150
• AA-19A-20CK	2 Fiber Optic Handpiece Connections (6 pin)\$1,100
• AA-20AMAL-CK	Amalgam Separator\$650
• AA-20A	Height Adjustable Pivoting Tray\$395
• 120367	Integrated AEU-5000 Electric Motor\$1,995

Recommended Accessories

• AHP-50	Lowspeed Air Motor (pg 22)\$325
• AHP-96	Highspeed Air Handpiece (pg 22)\$395
• AHP-101	1:1 Straight Handpiece (pg 22)\$275
• AHP-63MBFO-XL	1:1 Contra Angle (pg 11)\$1,195
• AHP-72MBFO-XL	1:5 Increaser Contra Angle (pg 11)\$1,295
• Scaler Tips	Hygiene and Periodontics (pg 13)\$175 ea

INTRODUCING THE AEU-7000L-70V ELITE . . . \$4,495**Highest-End Functionality at an Unbeatable Price!****Recommended
Handpiece****AHP-85MBFO-CX**
20:1 Fiber Optic
Contra Angle

- All the great features of the AEU-7000E-70V (see page 4), plus LED handpiece illumination*
- Compatible with most brands of fiber optic and non-fiber optic 1:5, 1:1, and 20:1 E-type contra angles.
- LED intensity is fully user-adjustable and savable to presets
- Virtually maintenance free! Never change a bulb or worry about delicate fiber optic cabling.

AE-230L-40 Replacement Micromotor \$875**

*Handpiece illumination requires AEU-7000L series motor and appropriate fiber optic handpiece

**AE-230L-40 micromotor is incompatible with AEU-6000 series, AEU-7000 & AEU-7000E models.



AEU-7000E-70V . . . \$3,995

Implant, Endo & Oral Surgery Motor - Top-End Functionality!

- Variable ratio selection: 1:5, 1:2, 1:1, 8:1 and 20:1 E-type handpieces as well as any reduction ratio that can be read from dynamometer calibration
- Advanced dynamometer calibration system accurately measures the properties of the micromotor and handpiece prior to use for ratios between:
 - Implant Mode 16:1 to 32:1
 - Endo Mode 4:1 to 16:1
- Six programmable preset buttons for storing and labeling custom settings in memory
- Upgradeable software for a longer return on investment
- Fully integrated, easy-load irrigation pump
- Adjustable torque:
 - Implant Mode will reach torque up to 60 Ncm, which enables a single, high-efficiency 20:1 reduction handpiece (AHP-85MB) to be used for the complete implant procedure
 - Endo Mode will reach torque up to 1000 g-cm using 8:1 reduction handpiece
- Switch between Endo Mode and Implant Mode at the touch of a button
- Powerful, brushless, 40,000 rpm, autoclavable micromotor
- Selectable Auto-Stop function in Implant Mode and Auto-Stop-Reverse in Endo Mode when desired torque setting is reached
- Contemporary design with easily cleaned surfaces
- Large, bright, easy-to-read display
- Intuitive user interface that features adjustable handpiece ratio, speed, torque, and irrigation flow settings
- Autoclavable motor holder may be free-standing or mounted to either side of the console
- Includes AE-70V Multi-Function Foot Control with Pump On/Off, Flow Rate, Micromotor Direction, Preset Selection, and Torque Cycle buttons
- Includes autoclavable tubing set for standard irrigation bags; optional tubing set accommodates rigid irrigation bottles
- Auto-sensing global voltage compatibility
- FDA, UL, and CE compliant
- Manufactured in the USA



Carrying Case

(For All AEU-7000 & AEU-6000 models)

410189 \$450.00

Replacement Parts

AE-70V-10 Foot Control with 10' Cable	\$500
AE-7PM (Optional) On/Off Footswitch	\$125
AE-23 Replacement Tubing Set	\$55
AE-23-BOT Tubing Set for Rigid Irrigation Bottles	\$55
AE-230M-40 Replacement Micromotor.	\$795
AE-230M-40EXT Micromotor with 10' Cable	\$895
461540 Dynamometer Adapter (AEU-7000 only)	\$8.50



AEU-6000-70V . . . \$2,850

Implant, Endo & Oral Surgery Motor - Exceptional Value!

- Ratio selection:
 - Implant Mode: 1:5, 1:2, 1:1, 20:1
 - Endo Mode: 8:1
- Calibration technology for accurate speed and torque performance
- Five programmable preset buttons for storing settings in memory
- Upgradeable software for a longer return on investment
- Fully integrated, easy-load irrigation pump
- Adjustable torque:
 - Implant Mode will reach torque up to 50 Ncm, which enables a single, high-efficiency 20:1 reduction handpiece (AHP-85MB) to be used for the complete implant procedure
- Endo Mode will reach torque up to 1000 g-cm using 8:1 reduction handpiece
- Automatically switches between Endo Mode and Implant Mode based on selected handpiece ratio
- Powerful, brushless, 40,000 rpm, autoclavable micromotor
- Selectable Auto-Stop function in Implant Mode and Auto-Stop-Reverse in Endo Mode when desired torque setting is reached
- Contemporary design with easily cleaned surfaces
- Large, bright, easy-to-read display
- Intuitive user interface that features adjustable handpiece ratio, speed, torque, and irrigation flow settings
- Autoclavable motor holder may be free-standing or mounted to either side of the console
- Includes AE-70V Multi-Function Foot Control with Pump On/Off, Flow Rate, Micromotor Direction, Preset Selection, and Torque Cycle buttons
- Includes autoclavable tubing set for standard irrigation bags; optional tubing set accommodates rigid irrigation bottles
- Auto-sensing global voltage compatibility
- FDA, UL, and CE compliant
- Manufactured in the USA

Recommended Handpieces

AHP-85MB-X 20:1 Reduction Contra Angle	\$1,395
AHP-85MB-CX 20:1 Reduction Contra Angle with Adjustable Depth Stop & Gauge	\$1,495
AHP-77D 1:2 Increaser Handpiece	\$1,195
AHP-71TI 1:5 Increaser Contra Angle	\$825
AHP-88MNP 8:1 Reduction Contra Angle Micro-Head for Endodontics	\$545

Did You Know?

The AEU-6000-70V meets or exceeds the functionality of other implant motors but sells for substantially less.

The motor cable and foot cable are also available in 10' lengths



(Handpieces sold separately)

AEU-260S . . . \$1,795

- Console displays handpiece ratios 1:5, 1:2, 1:1, 16:1 & 20:1
- Handpiece calibration
- Bright, easy-to-read display
- Intuitive user interface
- Selectable handpiece ratios
- Torque control (percentage)
- Auto-Stop
- Five programmable preset buttons
- Preset options allow speed, torque, direction and custom name preferences
- Powerful, brushless, 30,000 rpm, autoclavable micromotor
- AE-7PM 360° touch, weighted Footswitch
- Upgradeable software for a longer return on investment
- Auto-sensing global voltage compatibility
- Manufactured in USA

AE-70V Multi-Function Foot Control \$450

AE-225-30 Replacement Micromotor \$695

AEU-707Av2 \$1,995

- Easy-to-use
- Durable construction - thousands in use worldwide
- Powerful, brushless, 30,000 rpm, autoclavable micromotor
- Includes AE-20 Irrigation Tubing Set
- Console displays the following E-type handpiece ratios: 1:1, 16:1, 20:1, and 256:1
- Irrigation flow control, pump activation, handpiece ratio selection, micromotor direction, and speed adjustment
- Easily cleaned control panel
- Supplied with AE-7P On/Off Footswitch
- Dual voltage switchable - Specify 110 or 220 V when ordering
- Manufactured in the USA

Aseptico's Basic Implantology System with Irrigation

Excellent for General Dental Applications Too! Use with 1:5 Increaser Handpiece for Speeds Up to 150,000 RPM



(Handpieces sold separately)

AE-18A Variable Speed Foot Control Upgrade . . . \$500

- Hands-free pump activation, micromotor direction selection, and variable speed adjustment
- For use with AEU-707A or AEU-17B series motors
(No pump activation on AEU-17B series)

AE-4B-30 Replacement Micromotor \$495

Recommended Handpieces:

AHP-85MB-CX	20:1 Contra Angle \$1,495
AHP-85P-I	20:1 Contra Angle \$925
AHP-60NSK	256:1 Contra Angle \$1,095
AHP-77D	1:2 Increaser \$1,195
AHP-71TI	1:5 Increaser \$825

AE-7P On/Off Footswitch Replacement . . . \$125

Implant & Oral Surgery Systems Feature Comparison Chart

	LED Fiber Optics	Irrigation Pump	Brushless, Autoclavable Micromotor	Micromotor Speed	Dynamometer Handpiece Calibration	Automatic Handpiece Ratio Detection	Uses 1 Handpiece for Complete Implant Procedure	Digital Torque Control	Upgradeable Software	Programmable Presets	Displayed Handpiece Ratios	Variable Speed Foot Pedal	Endodontic Capabilities	Voltage
AEU-7000L-70V	✓	✓	✓	40,000 RPM	✓	✓	✓	✓	✓	✓	1:5, 1:2, 1:1 8:1, 20:1	✓	✓	Universal
AEU-7000E-70V		✓	✓	40,000 RPM	✓	✓	✓	✓	✓	✓	1:5, 1:2, 1:1 8:1, 20:1	✓	✓	Universal
AEU-6000-70V		✓	✓	40,000 RPM			✓	✓	✓	✓	1:5, 1:2, 1:1 8:1, 20:1	✓	✓	Universal
AEU-260S			✓	30,000 RPM			✓	✓	✓	✓	1:5, 1:2, 1:1 16:1, 20:1	OPTIONAL		Universal
AEU-707Av2		✓	✓	30,000 RPM							1:1, 16:1 20:1, 256:1	OPTIONAL		Dual Voltage Switchable



Handpieces for Implant Units

Model	Ratio	Max Torque	Speed Range <i>(Varies per machine)</i>	Price
AHP-85MB-X*	20:1	.80 Ncm	.15 to 2,000	\$1,395
AHP-85MB-CX*	20:1	.80 Ncm	.15 to 2,000	\$1,495
AHP-85MBFO-CX*	20:1	.80 Ncm	.15 to 2,000	\$1,895
AHP-85P-I	20:1	.55 Ncm	.15 to 2,000	\$895
AHP-85AS-I	20:1	.50 Ncm	.15 to 2,000	\$995
AHP-85-I	20:1	.50 Ncm	.15 to 2,000	\$895
AHP-63	1:1		.300 to 40,000	\$425
AHP-64	1:1		.300 to 40,000	\$500
AHP-77D (20° angle)	1:2		.600 to 80,000	\$1,195
AHP-77DS (straight)	1:2		.600 to 80,000	\$1,195
AHP-71TI	1:5		.1,500 to 200,000	\$825

Removable head for thorough cleaning
No tools required!



HANDPIECE INNOVATIONS!

*The Mont Blanc contra angles are made of stainless steel & titanium, making them very strong and lightweight. All Mont Blanc handpieces are easy to clean. They feature smooth surfaces for improved infection control, precision ball-bearing construction and their efficient output prevents overheating during most surgical procedures. The AHP-85MB features a fully detachable head, allowing for thorough maintenance to ensure lasting durability and reliability.

DTC® 2 High Performance Endodontic System Advanced Technology

Did You Know?

The AEU-26 presets can be customized for any rotary endo file system in the market.



(Handpieces
sold separately)

AEU-26 DTC® 2 . . . \$1,795

- Handpiece Ratios: 1:5, 1:1, 8:1, 16:1
- Advanced endodontic technology
- Handpiece calibration
- Bright, easy-to-read display
- Intuitive user interface
- Integrated file library
- Torque shown accurately in gram-centimeters
- Auto-Stop-Reverse reduces file separation
- Five programmable preset buttons
- Preset options allow single file, file series, or custom file series
- Auto-sensing global voltage compatibility
- Powerful, brushless, 30,000 rpm, autoclavable micromotor
- Upgradeable software for a longer return on investment
- Manufactured in USA



AE-70VE Multi-Function Foot Control . . . \$450

- On/Off Footswitch
- Buttons easily cycle through: Torque Modes, Presets, Motor Direction & Files in a Series

Replacement Parts

AE-225-30 Replacement Micromotor	\$695
AE-7PM On/Off Footswitch	\$125

Recommended Endo Handpieces

AHP-88MNP 8:1 Reduction Contra Angle Micro-Head for Endodontics	\$575
AHP-71TI 1:5 Increaser Contra Angle	\$825

Did You Know?

Aseptic also manufactures HandiDam®, our innovative, pre-framed rubber dam that is ideal for Endodontic procedures! See HandiDam and our other rubber dam products and accessories on pages 24-26.



AEU-25 Endo DTC® \$1,280

- Handpiece Ratios: 1:5, 1:1, 8:1, 16:1, 64:1
- Advanced digital technology that is easy to use
- Proven reliability
- Handpiece calibration for precision speed and torque control
- Torque shown accurately in gram-centimeters
- Auto-Stop-Reverse reduces file separation
- Five programmable preset buttons
- Integrated file library
- Powerful, brushless, 30,000 rpm, autoclavable micromotor
- Upgradeable software for a longer return on investment
- Includes AE-7P On/Off Footswitch
- Dual voltage switchable - Specify 110 or 220 V when ordering
- Use with peristaltic irrigation pump (AEU-7) for surgical applications
- Manufactured in the USA



AEU-7D Dual Foot Control \$180

Integrates on/off footswitch with second foot pedal to step through a file series or select preset with the AEU-25

Did You Know?

You can upgrade your AEU-25 for implant or general dental procedures just by adding the AEU-7 pump and an AHP-85P-I or AHP-71TI handpiece (see page 7)

Replacement Micromotors

AEU-25: AE-4B-30	\$495
AEU-20: AE-4B-10	\$495



AEU-20 Endo ITR® \$1,050

- Basic functionality
- Auto-Stop-Reverse reduces file separation
- Four levels of analog torque control
- Brushless, 10,000 rpm, autoclavable micromotor
- Console displays the following E-type handpiece ratios: 1:1, 8:1, 16:1, and 64:1
- Includes AE-7P On/Off Footswitch
- Specify 110 or 220 V when ordering
- Manufactured in the USA



AEU-7 Peristaltic Irrigation Pump \$775

- Superb electric-powered irrigation
- Includes AE-20 Autoclavable Tubing Set
- Compatible with AEU-25 Endo DTC® (pictured), AEU-20 Endo ITR, and AEU-17Bv2 SurgiMotor II
- Includes AE-7 Footswitch for independent operation
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AEU-7DC Foot Control \$180

- Simultaneous control of AEU-7 Irrigation Pump and either AEU-25 Endo DTC® or AEU-17B series of motors.

Handpieces for Endodontic Units

AHP-71TI	1:5 Increaser Contra Angle	\$825
	Achieves 150,000 rpm with 30K motor for canal access, friction grip	
AHP-63	1:1 Standard Contra Angle	\$425
Speed Range:	300 to 30,000 rpm with Endo DTC® 1,500 to 9,500 rpm with Endo ITR	
AHP-88MN	8:1 Reduction Contra Angle	\$495
	Mini-head for Endodontics	
AHP-88MNP	8:1 Reduction Contra Angle	\$545
	Micro-head for Endodontics	
AHP-62MN	16:1 Reduction Contra Angle	\$595
	Mini-head for Endodontics	
AHP-62MNP	16:1 Reduction Contra Angle	\$645
	Micro-head for Endodontics	
8:1 Speed Range:	40 to 3,600 rpm with Endo DTC® 100 to 950 rpm with Endo ITR	
16:1 Speed Range:	20 to 1,800 rpm with Endo DTC® 50 to 600 rpm with Endo ITR	



AHP-71TI



AHP-63



AHP-88MN
AHP-62MN



AHP-88MNP
AHP-62MNP



AHP-09 AseptiSpray™

A must for proper handpiece cleaning and maintenance! Lubricant and nozzle (included) are autoclavable.

AHP-09*	1 can (8.8oz)	\$19.95
	3 cans (26.4oz)	\$55.50
AHP-09N	Nozzle only	\$7.50

* Contents under pressure.

This product may only be shipped by surface.

Endodontic Systems - Feature Comparison

Features	AEU-26 DTC-2	AEU-25 Endo DTC®	AEU-20 Endo ITR
Brushless, Autoclavable Micromotor	✓	✓	✓
Motor Speed	30,000 rpm	30,000 rpm	10,000 rpm
Handpiece Calibration	✓	✓	
Auto-Stop-Reverse	✓	✓	✓
Internal File Library	✓	✓	
Customizable File Series	✓		
Digital Torque Control (g-cm)	✓	✓	
Analog Torque Control			✓
Upgradeable Software	✓	✓	
Programmable Presets	✓	✓	
Display Handpiece Ratios	1:5, 1:1, 8:1, 16:1	1:5, 1:1, 4:1, 8:1, 16:1, 64:1	1:1, 8:1, 16:1, 64:1
Dual Foot Pedal Option	Multi-Function Option	✓	
Voltage	Universal Voltage	Dual Voltage Switchable	110V or 220V

AEU-5000**Electric Operative Motor . . . \$1,995**

- Compact, easy-to-clean console
- Powerful, brushless, 40,000 rpm, autoclavable micromotor
- Innovative LED technology yields the brightest handpiece illumination in the industry!
- Fiber optic intensity adjustment
- Endo functionality
- Digital torque control
- Large, easy-to-read, back-lit display
- Intuitive user interface with large buttons
- Compatible with E-type handpieces for highspeed, lowspeed, Endo, and a variety of other dental applications
 - Increaser 1:5, 1:2
 - Standard 1:1
 - Reduction 5:1, 8:1
- Enable or disable any of the above handpiece ratios so that only user-desired handpieces are selectable on the display. No other electric operative motor allows this level of customization!
- Air/water irrigation through the micromotor and 1:5 increaser, 1:1 and 5:1 reduction handpieces
- Five programmable presets per handpiece ratio
- Attach to any Midwest 4-hole air connection; uses your disk-type foot control for operation
- Optional mounting bracket

AE-240-40 Replacement Micromotor \$695



330569 Adjustable Mounting Bracket . . . \$35

- Easily mount AEU-5000 console under a counter
- Viewing angle adjustable by 20°



AHP-72MBFO-XL



AHP-63MBFO-XL



AHP-56MBFO-XL

New Mont Blanc Handpieces

Model	Ratio	Internal Irrigation	Fiber Optics	Price
AHP-72MBFO-XL	1:5	✓	✓	\$1,295
AHP-72MB-X	1:5	✓		\$1,195
AHP-63MBFO-XL	1:1	✓	✓	\$1,195
AHP-63MB-X	1:1	✓		\$1,095
AHP-56MBFO-XL	5:1	✓	✓	\$1,195
AHP-56MB-X	5:1	✓		\$1,095

- Extremely lightweight due to integration of titanium
- Smallest head in the market
- Ceramic ball bearings for reduced heat - 1:5 only
- Special coating and new ergonomic design for improved handling
- Water flow of up to 80 ml/min
- Integrated spray allows efficient cooling of the bur, tooth and the contra angle itself, reducing risks of overheating and avoiding excessive wear and tear to the gear systems

**ALSO GREAT
FOR THE
MOBILE
HYGIENIST!**



Model ADU-17X2
(Handpieces sold separately)

ADU-17X Express II \$3,795

ADU-17X2 Express Plus II . . . \$3,995

Features:

- 3-way air/water syringe
- HVE and saliva ejector suction
- 1 Liter water reservoir
- Extra-large 1 gallon (3.79 liter) waste container with automatic overflow shut-off
- Built-in compressor
- Handpieces and suction may be activated via the attached footswitch or manually switched on to run continuously
- AWS-1C water connection for **ASC-12** or other external scalers
- All components pack into impact-resistant carrying case with wheels for easy transportation
- Convenient instrument holder attaches to the retractable handle on the case
- **Model ADU-17X** includes one highspeed pneumatic handpiece connection with water adjustment
- **Model ADU-17X2** includes two pneumatic handpiece connections (Highspeed has water adjustment)
- Dual voltage switchable - Specify 110 or 220 V when ordering - fuses for both voltages included
- Weight: 43.5 lbs (19.7 kg)
- Dimensions: 10.5" x 14.5" x 24" (27 x 36 x 61 cm)
- Manufactured in the USA

Options:

Handpiece Illumination System \$475
(See page 22 for details)

**Add ASC-12 Amplitude Piezoelectric
Scaler (see back cover) for only \$800
A \$195 Savings!**

(tips sold separately) **PROMO CODE 0211**

Recommended Handpieces:

- **AHP-96** Highspeed Air Handpiece (4-hole) \$395
- **AHP-50** Lowspeed Air Motor. \$325
- **AHP-101** 1:1 Straight Handpiece \$275




410136 ASC-12 Scaler Case . . . \$110

Add the ASC-12 Scaler (see back cover) to your portable equipment using this impact-resistant travel case.

AEU-425 Transport II \$5,495**AEU-425FO Transport II
with Fiber Optic Micromotor. . . \$5,995**

- Fully self-contained electric dental system
- Easy to set up and operate
- Powerful, brushless, 30,000 rpm, autoclavable micromotor
- Console displays the following E-type handpiece ratios: 1:5, 1:1, 4:1, and 16:1
- Fiber optic micromotor (AEU-425FO) features internal irrigation. Compatible with most brands of increaser contra angles with light
- 3-way air/water syringe
- 1 Liter water reservoir
- Automatic activation of high volume and saliva ejector suction
- Ample waste containers with automatic shut-off and fill indicator
- Built-in compressor
- Includes AE-7P On/Off Footswitch
- Air pressure gauge
- All components pack into impact-resistant carrying case with wheels for easy transportation
- Dual voltage (auto sensing) - Specify 110 or 220 V when ordering - fuses for both voltages included
- Power cords for North America (110 V) and Europe (220 V) included
- Dimensions: 10.5"x 15.4" x 24" (26.7 x 36.8 x 61 cm)
- Weight: 53.5 lbs (24.3 kg)
- Manufactured in the USA

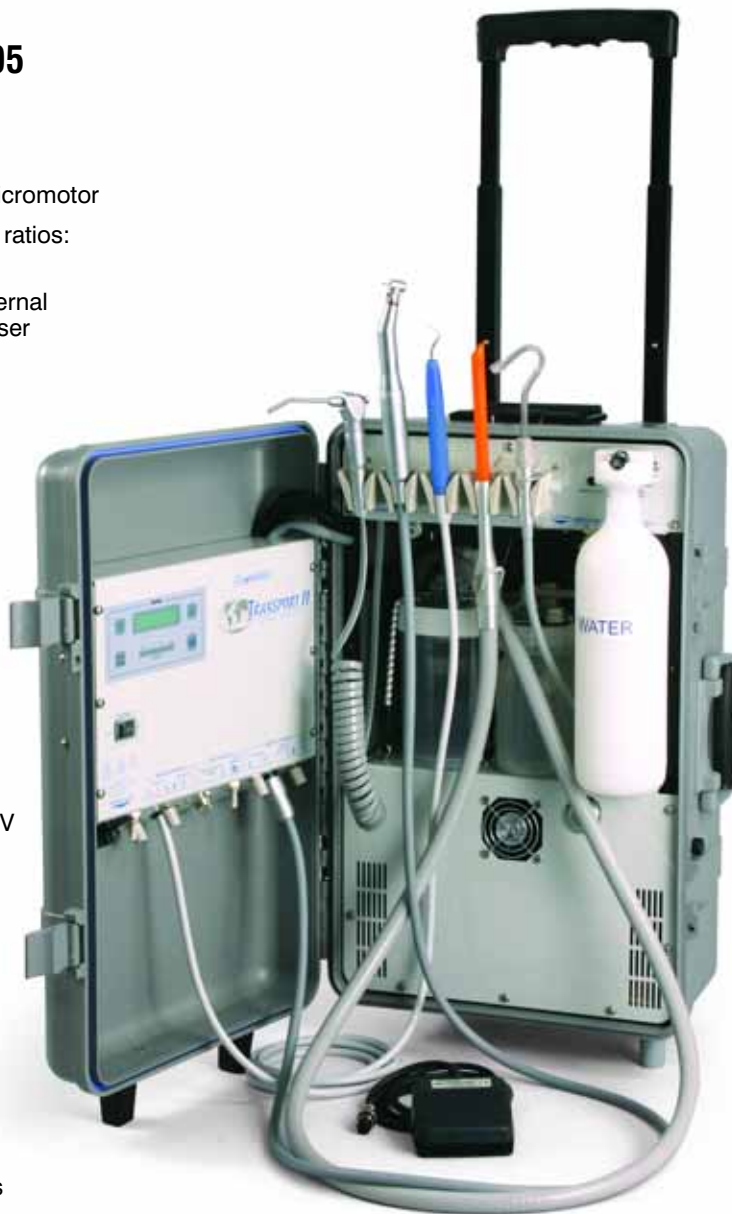**ASC-10-ACC Complete Kit . . . \$595**

includes 3 tips (*below*), 3 sleeves, 3 tip wrenches

ASC-10-PE37: Thin Sub Subgingival**ASC-10-PE38: Slim Universal****ASC-10-PE39: Power Universal**

Replacement Tips (ASC-10 only) \$175 each

Recommended Handpieces:

AHP-71TI (for AEU-425) Highspeed	\$825
AHP-72MBFO-XL (for AEU-425FO) Highspeed with Fiber Optics & Internal Irrigation	\$1,295
AHP-101 Lowspeed	\$275

Options:

ASC-10 Optional Piezo Ultrasonic Scaler (<i>pictured</i>)	\$1,095 (includes ASC-10-ACC complete kit)
AWS-1C Water Connection for ASC-12 (<i>see back cover</i>) or other external scalers.	\$75
AE-205-30 Replacement Micromotor for AEU-425	\$750
AE-200-30 Replacement Micromotor with Light for AEU-425FO	\$875
820039 Replacement Light Bulb for AE-200-30	\$60

ADU-12 Taskforce (Standard) . . . \$2,395

ADU-12DCE Taskforce Deluxe . . . \$2,750

- Fully self-contained economical dental system
- Easy to set up and operate
- 3-way air/water syringe
- Powerful saliva ejector suction
- Dual handpiece control with manual selection
- Compatible with Midwest 4-hole or Borden 2/3-hole pneumatic handpieces - specify when ordering
- Powerful internal oilless air compressor provides over 32 psi to handpieces and syringe
- 1 Liter water reservoir
- Ample 32oz waste container
- Includes On/Off Footswitch
- **ADU-12DCE Deluxe Model Includes:**
 - Handpiece off-switch to run suction independently
 - AWS-1C water connection for **ASC-12** or other scalers
- All components pack into sturdy carrying case for easy transportation
- Integrated combination lock for added security
- Dual voltage switchable - Specify 110 or 220 V when ordering - fuses for both voltages included
- Dimensions: 12.5"x 17.5"x 7.5"
(31.8 x 44.5 x 19.1 cm)
- Weight: 28.5 lbs (12.8 kg)
- Manufactured in the USA




**ADU-12 with AA-12A
& AA-12HVE options**

Recommended Handpieces:

- **AHP-55** Highspeed Air Turbine (Midwest 4-hole)
Maximizes air supply for superior cutting power . . \$295
- **AHP-50** Lowspeed Air Motor. \$325
- **AHP-101** 1:1 Straight handpiece \$275

**Add ASC-12 Amplitude Piezoelectric
Scaler (see back cover) for only \$800
A \$195 Savings!**

(tips sold separately) PROMO CODE 0211

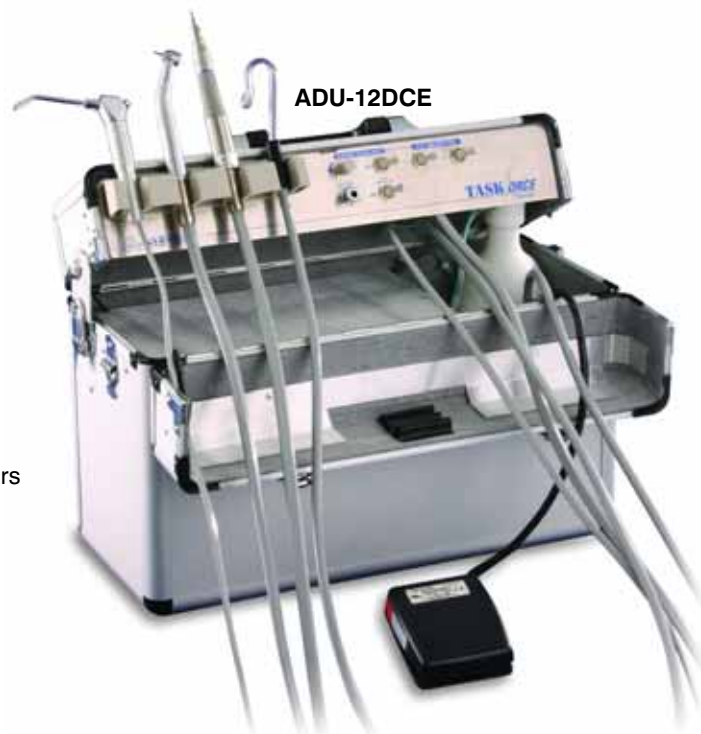

ADU-12DCE

Optional Accessories:

- **AA-12A** Handpiece Extension Arm - installed
(pictured with ADU-12) Conveniently positions
handpiece, syringe, and saliva ejector valve closer
to the operator \$55
- **AA-12HVE** High Volume Tip Adapter (pictured
with ADU-12) Enables use of standard 11mm
high volume tips for greater air flow \$15
- **AA-26SS** Stainless Steel Cuspidor with solids trap
5" (12.7mm) diameter for convenient expression
during treatment \$25
- Handpiece Illumination System. \$475
(See page 22 for details)



AA-12A



AA-26SS

Self Contained Dental Systems



Feature Comparison Chart

Features	AEU-425 AEU-425FO	ADU-17X ADU-17X2	ADU-12 ADU-12DCE
Electric Handpiece Operation	✓		
Pneumatic Handpiece Operation		✓	✓
High Volume Suction	✓	✓	
Saliva Ejector Suction	✓	✓	✓
Independent Handpiece & Vacuum Operation	✓	✓	✓ ADU-12DCE
Fiber Optics Standard	✓ AEU-425FO		
Optional Fiber Optics		✓	✓
Optional Ultrasonic Scaler	✓	✓	✓
Air/Water Syringe	✓	✓	✓
1 Liter Water Supply	✓	✓	✓
Water Quick-Connect for External Ultrasonic Scaler	OPTIONAL	✓	✓ ADU-12DCE
Optional HVE Tip Adapter			✓
Optional Handpiece Extension Arm			✓
Optional Stainless Steel Cupidor	✓	✓	✓
Built-in Combination Lock			✓
Padlock Loops (<i>padlock not included</i>)	✓	✓	
Regulatory Approval for Use in Europe and Canada	✓	✓	✓ ADU-12DCE Europe Only
Dimensions	10.5" x 14.5" x 24" 26.7 x 36.8 x 61cm	10.5" x 14.5" x 24" 26.7 x 36.8 x 61cm	12.5" x 17.5" x 7.5" 32 x 44.5 x 19cm
Weight	53.5 lbs (24.3 kg)	43.5 lbs (19.7 kg)	28.5 lbs (12.8 kg)
Voltage	Dual Voltage - Auto Sensing	Dual Voltage - Switchable	Dual Voltage - Switchable



Prestige Medical Classic Portable Sterilizer

ASU-12-110 (110V) \$1,245

ASU-12-230 (230V) \$1,395

- Compact and portable
- Easy to Operate
- 18 min sterilization cycle (23 min total cycle time)
- Takes only 6 minutes to reach sterilization temperature
- Temperature: 121°C (250°F)
- Lightweight: Only 11.5 pounds.
- Cost effective
- 9 liter capacity
- Includes instrument basket
- Depressurization valve to reduce cooling time.

The Only FDA, CE, & UL Certified X-rays for Handheld Use!



**ARU-07 NOMAD™ PRO
with Transport Case \$7,960**

**ARU-07NC NOMAD™ PRO
without Transport Case \$7,495**

- Same features as the ARU-06, plus:
- Lighter-weight - 5.5lbs (2.5 kg)
- Color LCD screen
- Simple touch controls



NOMAD™ Transport Case

410171 (specify ARU-06 or ARU-07) \$500

Replacement Parts

ARU-06 NOMAD™ CLASSIC

- CP-0002** Battery Charger
(specify 110 or 220 V) \$150
- CP-0001** Battery Pack \$130
- AP-0012** Backscatter Shield \$425

ARU-07 NOMAD™ PRO

- AP-0025** Battery Handle \$320
- AP-0035** Battery Charging Cradle \$114
(power cord sold separately)
- CP-0062** 120/230 V Charger Cord \$65
(plug adapter required for 230 V)

Comprehensive radiation safety data for NOMAD™ is available upon request.



**ARU-06CE NOMAD™ CLASSIC
with Transport Case \$6,995**

**ARU-06CENC NOMAD™ CLASSIC
without Transport Case \$6,495**

- Ideal for portable and forensic dentistry
- Hand-held unit is designed for operator safety and patient acceptance
- Backscatter guard and internal shielding protect user
- Easy to operate with adjustable exposure timer (.01 - .99 seconds)
- Battery-powered, cordless
- 60 kVp, 2.3 mA. Compatible with digital sensors or traditional X-ray film
- Weighs just over 8 lbs (4 kg)
- System includes: NOMAD™ Portable X-ray unit, backscatter shield, battery charger, and two batteries
- Specify 110 or 220 V battery charger
- Manufactured in the USA

**ARU-06BS Deluxe X-Ray Stand
for NOMAD™ X-Ray . . . \$795**

- Secure, telescoping frame
- Rotating X-ray mounting clamp
- Post extends from 36" to 60" (91.4 to 152.4cm)
- Includes X-ray remote switch
- Includes padded carrying case.
- Easy positioning of the X-ray with patient in chair is enhanced when the optional caster kit is installed (Not for use with Nomad™ PRO)



AA-01CK Optional Caster Kit \$225



Portable Lights

- Superior illumination in remote field environments
- Height-adjustable tripod stand with gooseneck
- Price includes padded carrying case
- Folded Dimensions: 9" x 45" x 5.5" (23 x 114 x 14 cm)
- Manufactured in the USA

410084 Replacement Padded Carrying Case\$125



**NEW
LED Light!**

ALU-40LED **\$975**

- Adjustable high Intensity LED light with seven beam diameters
- Replaceable protective lens cover
- Amber filter for working with composites
- Super long life & durable LED
- No bulbs to replace - EVER!
- 1,500 FC @ 24" (16,145 LUX @ 61 cm)
- Universal voltage input (100-240 VAC 50-60 Hz)
- Auxiliary DC input (12-24 VDC)
- Weight: 12 lbs (5.4 kg)

461936 Replacement protective lens cover. . \$12

AA-40DC-BC DC Battery-Power Clip Assy . . \$50

AA-40DC-LP Automobile Power Adapter . . . \$40

AA-40FK Color Filter 5-Pack \$10

AA-40LK Lens Replacement Kit..... \$120



ALU-27 **\$895**

- Adjustable beam diameter, light intensity adjustment, and power switch located at the lamp head
- Includes extra bulb
- 1,600 FC @ 24" (17,200 LUX @ 61 cm)
- Weight: 16 lbs (7.27 kg)
- Specify 110 or 220 V when ordering

820003 Replacement Bulb \$35

ALU-29 **\$875**

- Light intensity adjustment and power switch located on tripod stand
- Includes extra bulb
- Weight: 15 lbs (6.81 kg)
- 1,100 FC @ 22" (11,800 LUX @ 55.9 cm)
- Specify 110 or 220 V when ordering

820001 Replacement Bulb \$24

AA-98 Airpac II Compressor . . . \$1,695

- ¾ HP oilless compressor; adjustable to 100 PSI
- Compatible with all Aseptico pneumatic delivery systems
- Air pressure switch
- Sturdy aluminum chassis
- Built-in wheels and pull handle
- 1.0 gallon air tank with drain (3.79L)
- Quick connect for air line
- Displacement: 4.1 CFM (116 LM)



- Electrical: 10.5 Amps, 830 Watts, 110-120 V 60 Hz
- Air Delivery:
2.12 CFM at 100 PSI (60 LM at 689 KPa)
2.62 CFM at 50 PSI (74.2 LM at 345 KPa); at 60 Hz
- Dimensions: 8.5" x 12.75" x 21.5" (22 x 32 x 55 cm) (dimensions do not include handles)
- Weight: 41 lbs (18.5 kg)
- Manufactured in the USA

ADC-01 AseptiChair . . . \$1,495

Aseptico's World-Renowned Portable Patient Chair

- Easy-to-clean, durable upholstery (gray)
- Light-weight, only 30 lbs (13.5 kg)
- Sturdy, high-strength rolled steel construction
- Safe for patients up to 500 lbs (225 kg)
- Backrest adjusts from fully upright to supine
- Seat height is adjustable to five (5) positions, between 11" to 22" (28 to 56 cm)
- Toe-board folds down enabling the ADC-01 to be used as an examination chair
- No assembly required!
- Price includes padded case, arm slings, and headrest cushion – no hidden charges and no loose parts to misplace!
- Folded dimensions: 23" x 34" x 8" (58 x 86 x 20.3 cm)
- Manufactured in the USA

410193 Replacement Padded Carrying Case . . . \$195



ADC-08 AseptiStool \$595

- Made and upholstered with the same durable materials as the ADC-01 AseptiChair
- Suitable for doctors and assistants
- One-piece design sets-up or stores in seconds, not minutes!
- Price includes padded case – no hidden charges and no loose parts to misplace!
- Excellent mobility/swivels on four casters
- Adjustable height
- Folded dimensions: 18" x 22" x 8" (48 x 53 x 23 cm)
- Weight: 15 lbs (6.75 kg)
- Manufactured in the USA

410093 Replacement Padded Carrying Case . . . \$105



**ATC-09 Fold-A-Way Cart****\$225**

- Ideal for reducing back-strain while transporting equipment
- Quick and easy set-up
- Ample cargo tray means less trips to and from the work site
- Weighs only 10 lbs (4.5 kg)
- Folded dimensions 2.5" (6.35 cm) height; conveniently stores under car seat or in trunk
- Accommodates weights up to 275 lbs (125 kg) on flat surfaces; 110 lbs (50 kg) up stairs or curbs
- Suitable for transporting file boxes or supplies within traditional dental offices

ADU-04ST Portable Delivery System**\$1,350**

- Compatible with standard Midwest 4-hole or Borden 2/3-hole pneumatic handpieces. Specify connector type and tubing style (coiled or straight) when ordering
- 3-way air/water syringe
- One highspeed line with water flow control; one lowspeed line without water
- Quick connect to external air source, such as a compressor or air tank
- Dual handpiece control with manual selection
- 1 Liter water reservoir
- Disc-type foot control
- Air pressure gauge
- Easy to set up and operate
- Manufactured in the USA
- Includes stainless steel instrument tray
- Price includes padded carrying case and AA-96EXT 25 ft. (7.6 m) Extension Air Line
- Weight: 18 lbs (8.1 kg)

Optional Accessories:

- **AA-21ST** High Volume Vacuum (*pictured*) . . . \$475
- **AA-25ST** Low Volume Vacuum \$400
- **AWS-1C** Water Connection for ASC-12 (*see back cover*) or other external scalers . . . \$75
- Handpiece Illumination System (*installed*) (*see page 22*) \$475
- Air-drive Handpieces (*see page 22*)

*(Includes Carrying Case)***AA-20 Portable Saliva Ejector . . . \$645**

- Powered via rechargeable battery pack, wall outlet or 12 VDC vehicle source
- Recharge battery via 110 to 220 V wall outlet or 12 VDC car lighter power adapter (*included*)
- Adjustable 80 - 550 mm Hg vacuum range
- 800 cc waste canister with float shut-off
- Dimensions: 9" x 7" x 8" (22.9 x 17.8 x 20.3 cm)
- Weight: 6.3 lbs (2.84 kg)
- Manufactured in the USA

AA-20BF - Replacement Bacteria Filters

Pkg. of 12 \$195

AA-20CAN - Replacement 800ml canister \$25**Optional Accessories:**

- AA-12HVE \$15
- AA-26SS cuspidor \$25

**ADU-01C Supply Case****\$195**

- Sturdy EMT carrying case for dental/medical supplies
- Removable, adjustable drawers
- Padlock loops (*padlock not included*)
- Dimensions: 10.75" x 19.5" x 15" (27 x 50 x 38 cm)
- Weight: 18 lbs (8.1 kg)

ADU-20B Dental Cart System . . \$2,595

- Improved chassis design
- Automatic dual handpiece selection with individual air and water adjustments
- Water on/off and flow adjustment for highspeed handpiece
- 3-way air/water syringe
- 1 Liter water reservoir
- Wet-dry foot control
- High volume vacuum and saliva ejector suction
- Compatible with Midwest 4-hole or Borden 2/3-hole pneumatic handpieces. Specify connector type and tubing style (coiled or straight) when ordering
- Air pressure gauge
- Height adjustable: 23" to 32" (60 to 81 cm)
- U-shaped base with casters for easy mobility between rooms
- Includes Ritter-size stainless steel tray: 13.5" x 9.8" (34 x 25 cm)
- Quick connect to air source – no junction box required!
- Compatible with Handpiece Illumination System (see page 22)
- Manufactured in the USA

Optional:

AWS-1C Water Connection for ASC-12 (see back cover)
or other external scalers \$75

**Add ASC-12 Amplitude Piezoelectric
Scaler (see back cover) for only \$800
A \$195 Savings!**

(tips sold separately) PROMO CODE 0211

**ADU-22 Dental Cart System . . . \$1,695**

- Dual handpiece control with manual selection
- 3-way air/water syringe
- Water on/off and flow adjustment for highspeed handpiece
- 1 Liter water reservoir
- Disc-type foot control
- High volume vacuum
- Compatible with Midwest 4-hole or Borden 2/3-hole pneumatic handpieces. Specify connector type and tubing style (coiled or straight) when ordering
- Height adjustable: 23" to 38" (59 to 97 cm)
- Star-shaped base with casters for easy mobility between rooms
- Includes Ritter-size stainless steel tray: 13.75" x 9.875" (35 x 25 cm)
- Quick connect to air source – no junction box required!
- Compatible with Handpiece Illumination System (see page 22)
- Manufactured in the USA



ATC-06 Tray Cart Stand \$620**ATC-04 Cart Stand \$275**

- Keeps Aseptico Endo or Implant motors organized and resolves operatory space concerns
- Large, removable Mayo tray on 24" (61 cm) swing arm (ATC-06 only)
- Height adjustable: 23" to 38" (59 to 97 cm)
- Star base with casters allows for easy mobility between rooms

(Surgical systems sold separately)

**ATC-08A Ritter Tray Cart \$395****ATC-08 Mountable Ritter Tray Cart . . . \$415**

- Large, removable Ritter-size tray
- Height adjustable: 23" to 38" (59 to 97 cm)
- Star base with casters allows for easy mobility between rooms
- Manufactured in the USA

Recommended Options for ATC-08

ADU-07A Triple Handpiece Control.	\$995
(pictured)	
ADU-07 Dual Handpiece Control	\$695

**ATC-02 Double Mayo Tray Cart . . \$695****ATC-01 Mayo Tray Cart \$425**

- Large, removable Mayo tray on top
- ATC-02 has additional tray on a 24" (61 cm) swing arm
- Height adjustable: 23" to 38" (59 to 97 cm)
- Star base with casters allows for easy mobility between rooms

**AA-50 Ritter size Tray . . \$55**

13.75" x 9.875" x 0.75"
(35 x 25 x 2 cm)

AA-51 Mayo Tray \$85

19.75" x 12.63" x 0.75"
(48 x 33 x 2 cm)



Air-Drive Handpieces

- **AHP-96** Highspeed air handpiece (4-hole) \$395



- **AHP-97** Highspeed air handpiece (5-hole) with fiber optics \$595



- **AHP-55** Highspeed air handpiece (4-hole)
Recommended for Taskforce series \$295



- **AHP-50** Lowspeed air motor (4-hole)
Adapts air turbine to accept E-style handpieces \$325



- **AHP-101** 1:1 Straight handpiece with post for AHP-53 Doriot-style handpiece. \$275



- **AHP-53** Doriot-style. Accepts latch-type burs. Connects to AHP-101 \$95

AHP-53 & AHP-101 Configuration**Did You Know?**

The AHP-63 (see page 10) can be used instead of the AHP-53 w/ AHP-101

USA customers can make all of your purchases on our Online Store!

www.aseptico.com

FiberBrite Handpiece Illumination System

AA-19A-04 (5-hole) \$475

AA-19A-06 (6-pin) \$475

A compact handpiece illumination system that installs easily to any handpiece control and is totally compatible with any standard fiber optic handpiece. An easy to replace halogen lamp provides brilliant, cool light for intra-oral use.

- ISO standard fiber optic handpiece terminal with stainless steel adapter and nut (Midwest 5-hole)
- Halogen lamp, 3.5 VAC, 720 mA
- Power transformer, 120 VAC, output 3.5 VAC, 700 mA
- Air electric switch activated
- Available in straight gray tubing only
- Specify 110 or 220 V when ordering
- Includes installation in Aseptico units

AA-19AB Extra Bulb Module (5-hole) \$65



AA-19A-04



ADU-04 Two-Handpiece Manual Control System \$695

- Dual handpiece control with manual selection
- Ideal for arm or wall mounting
- 3-way air/water syringe
- Water on/off and flow adjustment for highspeed handpiece
- Disc-type foot control
- Requires connection to air and water sources
- Compatible with Midwest 4-hole or Borden 2/3-hole pneumatic handpieces. Specify connector type and tubing style (coiled or straight) when ordering
- Manufactured in the USA



ADU-07 Economy Air System . . . \$695

ADU-07A Economy Air System . . \$995

- Flat-top design for mounting into cabinets
- **ADU-07A:** Triple handpiece control with automatic selection and individual air and water adjustments
- **ADU-07:** Dual handpiece control with manual selection
- Air pressure gauge
- 3-way air/water syringe
- Water on/off and flow adjustment for highspeed handpiece
- Disc-type foot control
- Requires connection to air and water sources
- Compatible with Midwest 4-hole or Borden 2/3-hole pneumatic handpieces. Specify connector type and tubing style (coiled or straight) when ordering
- Manufactured in the USA

Clean Water Systems

NWS-7 \$195

- 1 Liter white plastic bottle (not autoclavable)
- Includes mini air pressure regulator, mounting hardware, and interconnect tubing
- Specify wall or 2" (5 cm) post-mount configuration

NWS-8 \$195

- 1 Liter white plastic bottle (not autoclavable)
- Designed for Aseptico carts and portable dental units
- Includes interconnect tubing



AA-05A Folding Arm \$295

- Folding arm for wall mounting Aseptico Economy Air Units
- Dimensions: 24" (61 cm) extension

SUPERIOR ELASTICITY & TEAR-RESISTANCE!



HandiDam®

The Original Pre-Framed Rubber Dam!

- NEW generation dental dam
- Meets FDA "Powder Free" standard
- Favored among key opinion leaders for quality and reliability
- Pre-framed, flexible design facilitates access to the oral cavity for suction, x-ray films, or digital x-ray sensors
- Extremely low protein content reduces patient irritation (<50 micrograms)
- Isolates single tooth or group of teeth for endo and other dental treatments
- Saves time - eliminates the need to remove and replace traditional dam during the procedure!
- Consistent uniformity throughout
- Ergonomic and comfortable for greater patient acceptance
- Available in Latex (blue medium) or Latex-Free (mauve heavy). Specify when ordering:
- Latex-Free now made of Isoprene - very similar performance to regular HandiDam



Regular:

I-70BM Pkg. of 20

\$27.95

I-70BM-100 Pkg. of 100

(Best value!)

\$115

Latex-Free:

(Now Made of Isoprene)

I-70LF Pkg. of 20

\$39.95

I-70LF-100 Pkg. of 100

\$185

I-70K HandiDam Kit . . \$375

Includes all of the essentials for the HandiDam technique:

- **I-19** Punch
- **I-14** Forceps (University of Washington style)
- **KB-08** Stainless steel clampboard (see page 25)
- 8 of Aseptico's most popular clamps (see page 26)
- 2 boxes of regular HandiDam (40 pieces total)



NEW GENERATION DENTAL DAM!

- Superior elasticity and tear-resistance
- Meets FDA "Powder Free" standard
- Prepackaged in sealed bags of 36 for longer shelf life
- Made from the **highest quality** pure, natural latex
- Extremely low protein content reduces patient irritation (<50 micrograms)
- Pleasant vanilla scent
- Consistent uniformity throughout
- Available in Latex or Latex-Free
- Latex-Free now made of Isoprene - very similar performance to latex dam
- 6" x 6" (15 x 15 cm) squares



Description	Box of 36	Box of 288
Blue, thin	I-65BT . . . \$10.95	I-65BT-288 . . . \$60
Blue, med	I-65BM . . . \$10.95	I-65BM-288 . . . \$65
Blue, heavy	I-65BH . . . \$11.75	I-65BH-288 . . . \$70
Green, thin	I-65GT . . . \$10.95	I-65GT-288 . . . \$60
Green, med	I-65GM . . . \$10.95	I-65GM-288 . . . \$65
Green, heavy	I-65GH . . . \$11.75	I-65GH-288 . . . \$70



Clamp Boards

Organize your dam clamps with these autoclavable holders.

- KB-08** - Holds 8 clamps. Size 3.5" x 3.5"
(8.89 x 8.89 cm) \$45
- KB-12** - Holds 12 clamps. Size 3.5" x 6"
(8.89 x 15.24 cm) \$55
- KB-18** - Holds 18 clamps. Size 4.75" x 6.25"
(12.07 x 15.88 cm) \$65

Dam Punch



- Punch (Ainsworth Style)
I-18 \$89.95



- Rubber Dam punch that is especially good for HandiDam use.
I-19 \$129

Dam Frames



- U-shaped
- Lightweight
- Autoclavable stainless steel
- Two sizes available:
I-16 adult \$21.50
I-17 child \$18.50


Dam Forceps



- Forceps - Ivory Style
I-13 \$75



- Forceps - University of Washington Style
I-14 \$75

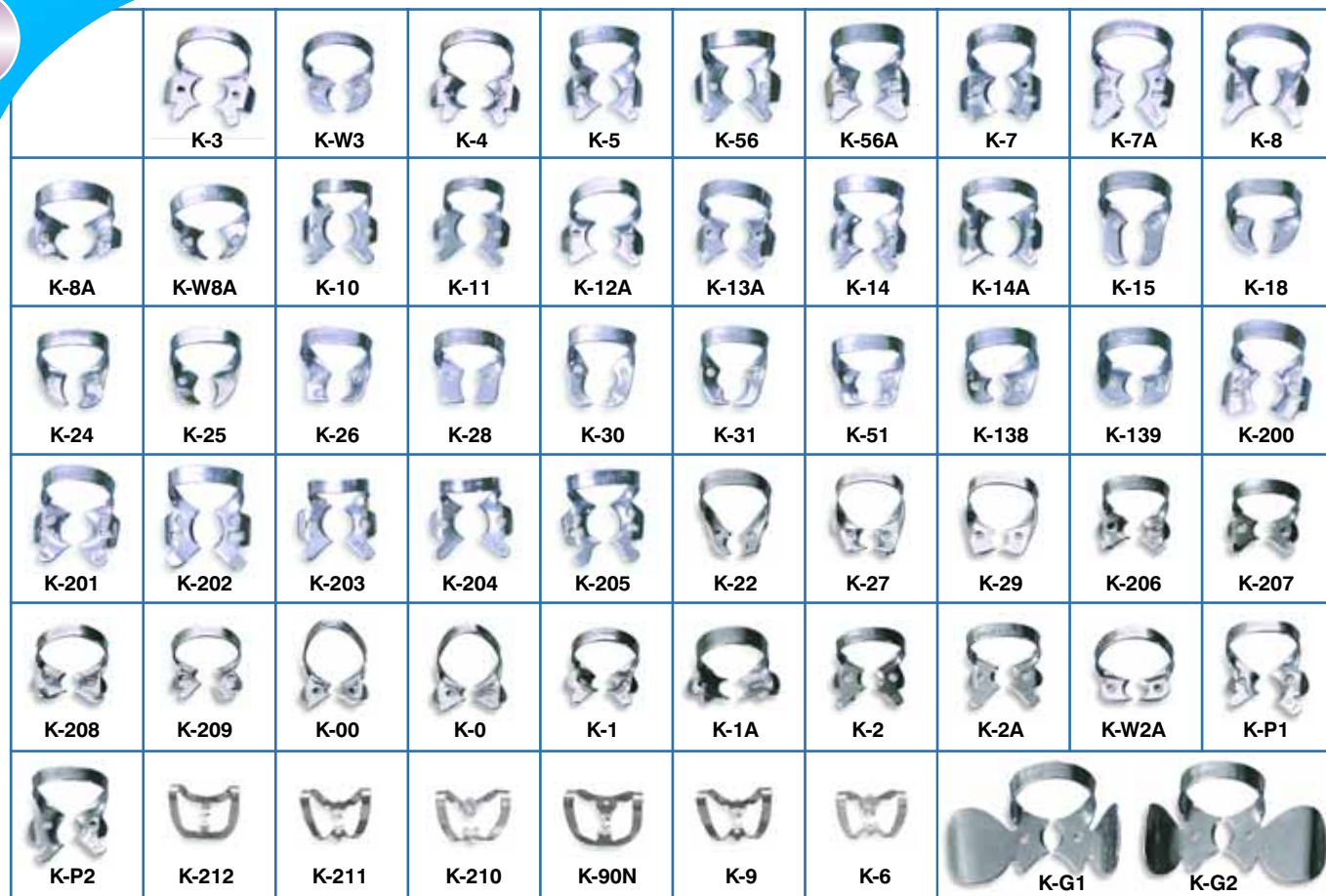

Rubber Dam & Labial Clamps \$10.75 each or 10 for \$100 • Gauze Clamps \$15.95 each

K-3For small lower molars, with flat beaks.
K-W3Similar to #3, but wingless.
K-4Similar to #3, but with festooned beaks.
K-5For lower molars.
K-56Similar to #5, but with festooned beaks.
K-56AFor larger lower molars, similar to #56 but with larger beaks.
K-7For lower molars, with flat beaks.
K-7AFor very big lower molars.
K-8For lower molars, with round beaks.
K-8AFor lower molars, with strong festooned beaks.
K-W8ASimilar to 8A, but wingless
K-10For small molars with short crowns - Left
K-11For small molars with short crowns - Right
K-12AFor third molars with serrated beaks - Left
K-13AFor third molars with serrated beaks - Right
K-14For half erupted molars.
K-14ASimilar to 14, but with large beaks.
K-15For big lower molars, with wide beaks, strong spring
K-18For upper and lower molars with round type beak.
K-24For deep buccal cavities of molars, wider fringe on buccal side - Left
K-25For deep buccal cavities of molar, wider fringe on buccal side - Right
K-26For upper molars, especially for erupted tooth
K-28For upper & lower molars, with very wide fringe for broad sight.
K-30For beaked molars - Left
K-31For beaked molars - Right
K-51For beaked molars, with both sides beaked
K-138For wisdom teeth, with stiff spring, serrated beaks - Left
K-139For wisdom teeth, with stiff spring, serrated beaks - Right
K-200For lower molars, with wings, broad fringe.
K-201Similar to #200, but with festooned beaks.

K-202For large molar clamps, similar to #200 but with large beaks
K-203For small lower molars & children's first molars - Left
K-204For small lower molars & children's first molars - Right
K-205For upper molars, with stiff spring.
K-22Similar to #207, but wingless
K-27Similar to #206, but wingless, festooned.
K-29For upper & lower bicuspid, with broad beaks for sider sight.
K-206For upper & lower bicuspid, with festooned beaks.
K-207For upper & lower bicuspid, with flat beaks.
K-208For big bicuspid, with similar pattern to #207
K-209For lower bicuspid, with flat beaks.
K-00For very small bicuspid & pedo centrals.
K-0For small bicuspid & pedo centrals
K-1For roots, with deep festooned beaks.
K-1AFor cervical cavities of bicuspid.
K-2For lower bicuspid, with flat beaks.
K-2ASimilar to #2, but with large beaks.
K-W2ASimilar to #2A, but wingless.
K-P1,P2For children's first molars.

LABIAL CLAMPS

K-210For centrals.
K-211Universal for labial cavities on anteriors.
K-212Cervical clamp for labial cavities, Ferrier.
K-6Similar to #9.
K-9Most popular anterior clamps.
K-90NStrong anterior clamps.

GAUZE CLAMPS

- Not for Rubber Dam
- Used with gauze or cotton rolls
- Extra-large wing acts as tongue press

**MOST POPULAR
CLAMPS IN BLUE**



Technique Chain & Grid

Technique Chain & Grid Traction easily attaches to a pre-emerged tooth, so that an orthodontic device may be used to draw the tooth into proper position. Aseptico is the originator of this remarkable traction system that professionals use worldwide!

Manufactured in the USA


AMA-09B Instrument . . . \$49

- Squeeze to release instrument allows for precise placement of grid
- Autoclavable



Chain & Grid Refill Kits . . . \$265

- 10 sets of Chain/Grid per kit.

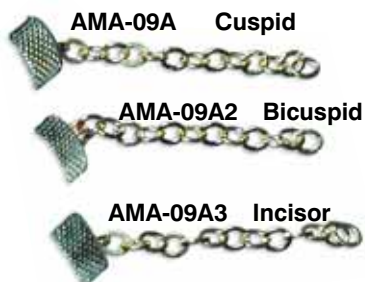

AMA-09D Lingual Buttons . . . \$35

- Six per package



AMA-09C Bonding Agent . . . \$105

- Non-etch glass ionomer bonding agent.

**AEU-16C Torque Plus \$795**

- Ideal for heavy-duty laboratory applications
- Powerful 35,000 rpm combined micromotor/handpiece
- Variable speed adjustment
- Forward/Reverse switch
- Selectable manual or foot control operation
- Includes AE-7P On/Off Footswitch
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AE-5L-32 Replacement Micromotor/Handpiece \$350**AEU-12C Torque Plus \$995**

- Same great features as the AEU-16C console!
- Ideal for medium or heavy-duty laboratory applications
- Powerful 32,000 rpm micromotor accepts E-type handpieces
- Includes AHP-64 straight 1:1 handpiece
- Includes AE-7P On/Off Footswitch
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AE-5C-30 Replacement Micromotor \$250**AHP-64 Replacement 1:1 Handpiece \$500****AEU-08D Laboratory Motor \$695**

- Ideal for medium-duty laboratory applications
- Powerful 30,000 rpm micromotor accepts E-type handpieces
- Includes AHP-101 1:1 E-type Handpiece
- Variable speed adjustment
- Forward/Reverse switch
- Selectable manual or foot control operation
- Includes AE-7P On/Off Footswitch
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AE-5C-30 Replacement Micromotor \$250**AHP-101 Replacement 1:1 Handpiece. \$275****AEU-10SS Power-Tip \$595**

- Ideal for light-duty laboratory applications
- 20,000 rpm micromotor accepts E-type handpieces
- Includes AHP-101 1:1 E-type Handpiece
- Variable speed adjustment
- Forward/Reverse switch
- Selectable manual or foot control operation
- Includes AE-7P On/Off Footswitch
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AM-103A Replacement Micromotor \$149**AHP-101 Replacement 1:1 Handpiece. \$275****AE-14A Variable Speed Foot Control Upgrade \$195**

- Compatible with AEU-10SS, AEU-12C, AEU-16C

AEU-03SS Porta-Tip \$695

- Ideal for light-duty laboratory applications
- Battery-powered for greater versatility
- 16,000 rpm micromotor accepts
- Includes AHP-101 1:1 E-type Handpiece
- Variable speed adjustment
- Specify 110 or 220 V when ordering
- Manufactured in the USA

AM-103A Replacement Micromotor \$149**AHP-101 Replacement 1:1 Handpiece. \$275**



BA-8 Autoclavable Bur Blocks

- Hold FG, HP, or LA burs in colored, tapered inserts
- Blocks can be connected together.
- Patented design
- Specify color when ordering: blue, yellow, white, orange, green, or red

Quantities: 1-5 \$8.95 ea.
 6-11 \$8.45 ea.
 12 or more \$7.95 ea.



AI-50 Crown Placer/Remover . . . \$95

Stainless steel crown placer/remover with contoured rotating grippers that lock onto the crown, and "bulldog" grip hemostat locking ratchet, make the removal or placement of crowns easy. Measures approximately 6.5" and is fully autoclavable.

AI-51 Replacement autoclavable rubber cushions (pair) . . . \$60

Loupe Specifications

Diopter	Magnification	Working Distance
3	1.75X	12" - 14" (30.5 - 35.5 cm)
5	2.25X	8" x 10" (20.3 - 25.4 cm)
7	2.75X	5" - 7" (12.7 - 17.8 cm)

Clip-On Loupes

\$34.95

- Flip-up design with flexible wire bracket fits any glasses
- Lightweight with two pivot points for adjustment
- Durable coating resists scratches
- Specify diopter size when ordering:
G-3 (3 diopter) • **G-5** (5 diopter) • **G-7** (7 diopter)



Clip-On Loupes

\$44.95

- All the same great features of above loupes, but with alligator clips.
- Convenient for mounting the loupes on any eyeglasses.
- Specify diopter size when ordering:
G-3C (3 diopter) • **G-5C** (5 diopter) • **G-7C** (7 diopter)



Headband Loupes

\$39.95

- Lightweight, comfortable headband-mounted loupes
- Convenient flip-up design
- Does not obstruct view
- Specify diopter size when ordering:
G-1 (3 diopter) • **G-2** (5 diopter)



Visor Loupes

\$44.95

- Ideal for laboratory use
- Convenient flip-up design
- Specify diopter:
G-13 (3 diopter) • **G-15** (5 diopter) • **G-17** (7 diopter)



**Handpiece Terminals**

- 3-hole, Borden style
AA-91 \$20.00
- 4-hole, Midwest style
AA-90 \$20.00
- Midwest to Borden connection adapter
730166 (not pictured) \$75.00

**Autoclavable Vacuum Lever Valve - Universal**

- With easy quick connect and lever on/off. Fits A-Dec™ holders
AA-35LAPW \$65.00

**Autoclavable Vacuum Lever Valve - Standard**

- With easy quick connect and lever on/off. Fits AA-69 holder
AA-36LA \$60.00

**Autoclavable Saliva Ejector Lever Valve**

- With easy quick connect and lever on/off. Fits AA-68 holder
AA-37LA \$40.00

**Autoclavable Vacuum Valve O-ring Service Kits**

Serves one AA-36LA or AA-35LA valve. Includes 3 autoclavable and 1 non-autoclavable o-ring

- AA-35LR** \$5.00



Serves one AA-37LA valve. Includes 2 autoclavable and 1 non-autoclavable o-ring

- AA-37LR** \$4.25

**Autoclavable Vacuum Valve Replacement Levers**

For AA-35LA & AA-36LA HVE valves

- AA-35HA** \$5.50

For AA-37LA saliva ejector valves.

- AA-37HA** \$5.50

AE-60 Electric to Pneumatic Footswitch Conversion Kit

Activates AEU-17B, AEU-20, AEU-25, or AEU-707A series motors via a standard pneumatic foot control \$105.00

**Replacement Tip - Saliva Ejector**

- AA-37RTG grey** \$5.50

**Saliva Ejector Screens**

Wire Autoclavable Screens

- AA-38W pkg. of 2** \$20.00

Disposable Screens

- AA-38 pkg. of 25** \$10.75

**TA-90D 3-Way Air/Water with Quick Change Tip**

- TA-90D** \$67.50

- TA-90DVC** \$10.75
 Replacement valve cartridge

- TA-90DVK** \$7.95
 Valve repair kit

**Autoclavable Syringe Tips**

Ensure OSHA compliance - autoclave syringe tips. Tips fit most recent model syringes.

- TA-1** \$7.50 ea.
 6 for \$40.00
 12 for \$75.00

Specify quantity when ordering:
 1, 6, or 12

**Quick Disconnects**

1/4" I.D. female x 1/8" MPT shutoff

- AA-61** \$19.95

1/4" I.D. female x 1/4" MPT shutoff

- AA-65** (not pictured) \$19.95

3/8" I.D. female x 1/4" MPT shutoff

- AA-63** (not pictured) \$29.95

1/4" I.D. male x 1/4" poly, without shutoff

- AA-60** \$19.95

3/8" I.D. male x 3/8" poly, without shutoff, used for Stryker drills

- AA-62** \$25.00

Schrader medical gas male swivel connector. Connects to 1/4" MPT. Type used in most hospitals for Hall drills.

- AA-64** \$95.00

**Saliva Ejector Hose Assembly:**

Saliva Ejector Adapter: Adapts 3/8" O.D. saliva ejector tubing to 7/16" O.D. HVE tip.

- AA-37A** \$7.50



Saliva Ejector Hose Assembly: Adapts to 7/16" O.D. HVE tips. Includes saliva ejector tip, 3 ft. of tubing and adapter.

- AA-37R** \$21.75



Saliva Ejector Hose Assembly with autoclavable AA-37LA lever control valve. Adapts to 7/16" O.D. HVE tip

- AA-37X** \$60.00



AE-20

Tubing & Accessories

- Fully autoclavable silicone
- Sets include bayonet & inline connectors

AE-23 - For AEU-7000 & AEU-6000 series \$55.00

AE-23-PUMP - For AEU-7000 & AEU-6000 series . \$14.95

AE-23-BOT - For AEU-7000 & AEU-6000 series rigid irrigation bottles. \$55.00

AE-20 - For AEU-707A series surgical systems & AEU-7 pump \$55.00

AE-26 - For AEU-925 VCT \$55.00

AE-2 - "Y" Connector \$3.25

AHP-03 - Irrigation Clips (10). \$15.00

AE-3 - Inline Barb Connector \$1.50

6" (15.2 cm) Straight Tubing:

AHP-07T \$4.75 ea.

AE-10 - Irrigation Tubing \$2.75/ft.



AHP-07 Cannula Clip Set \$16.95

Specially designed set includes silicone tubing, handpiece body clip, contra angle clip, cannula and nipple (as pictured) for internal irrigating Implant burs. Clip set fits most E-type contra angles. Fully autoclavable.

AHP-04: Cannula & Tubing \$7.50

AHP-05: Small Clip \$5.50

AHP-06: Large Clip \$5.50

Handpiece Accessories

AHP-85MB-CAN - Mont Blanc Cannula \$29.95

AHP-85MB-CLP - Mont Blanc Irrigation Clip \$19.95

AHP-85MB-STP - Mont Blanc Stop Only \$143.75

AE-2MB - Mont Blanc Y Connector \$39.95

AE-8D - Water Clip for AHP-64. \$75.00

700113 - Water Clip for AHP-77D or AHP-77DS . . \$75.00

Ordering Policy:

- Orders under \$50.00 add \$5.00 handling fee.
- Normal payment is via credit card. Aseptico accepts Visa, MasterCard, Discover, and American Express.
- No split payments on credit card. Entire amount of order is taken at time of original billing.
- Open account is available only if credit is approved in advance.
- Freight costs are added to all invoices
- Delinquent accounts are subject to 1.5% monthly interest, attorney's fees, collection and court costs when needed to satisfy debt.
- Washington and Florida state residents – add sales tax based on destination address.
- Residents of Dallas & San Antonio Texas - add state and local sales tax; all other Texas locations add state tax only.
- Purchasers outside of Washington State are responsible for any applicable state or use taxes.
- Prices subject to change without notice.
- Service and Repair charges are billed on credit card. Customer is expected to accept charges up to \$300 without advance notification. Minimum service/repair charge is \$95.
- Please visit our website: www.aseptico.com/policy.html for additional warranty and repair policy information.

Return Policy:

- Product that has been in clinical use cannot be returned for credit.
- All returns require prior authorization. Contact Aseptico Customer Service for a Return Authorization number (RA#).
- The RA# must appear clearly on the exterior of the package being returned.
- Packages not so marked may be rejected and returned to the sender.
- Packages received without prior authorization are not the responsibility of Aseptico. Aseptico reserves the right to dispose of such packages and their contents 30 days after receipt with no obligation to the sender.
- The sender is responsible for proper packaging and for all costs for return shipment to Aseptico.
- In general, new product may be returned for full credit (less shipping charges) under a Return Authorization number within 30 days of date of original invoice.
- In exceptional situations and at Aseptico's discretion, certain items may be returned for credit after the 30 day period. A 25% restocking fee will apply.

Aseptico's Quality Management System is certified for full quality assurance by the Conformity Assessment Body Laboratoire national d'essais / G-MED under the Medical Device Directive 93/42/EEC Annex II.3 and also according to ISO 13485 under the Canadian Medical Device Conformity Assessment System (CMDCAS).



Made In The USA

Important Notes to International Customers:

Most Aseptico products are CE-marked and Health Canada licensed. However, customers should verify the status of any specific item(s) with their sales representative or local distributor prior to finalizing the purchase.

Please contact your local Aseptico distributor for pricing and availability. Prices shown in the catalog are in US dollars.

If you are unsure of the Aseptico distributor in your country, for specific ordering and payment policies, or for general product inquiries, please contact our International Sales Department at Int@Aseptico.com.



ALU-29CF Portable Light

- Bright halogen light on highly adjustable tripod stand
- 10 degree beam, long life lamp provides 1,100 ft. candles at working distance of 24" (56 cm)
- Dual voltage: 115/220 V, 50/60 Hz
- Includes soft padded carrying case
- Fits in ADU-01CS mil-spec. hard chair case
- Weight: 15 lbs (6.75 kg)
- Manufactured in the USA

NSN: 6520-01-446-4170



ADC-08CF COMMAND AIR Portable Operator Stools

- Interchangeable stools for doctor and assistant
- Comfortable padded seat and backrest
- Castored legs with adjustable height
- Two ADC-08 stools in soft carrying cases store in one hard mil-spec. case
- Very easy to fold up and carry
- Dimensions: 25.5" x 24" x 19", 6.73 cu. ft. (65 x 61 x 48 cm)
- Weight:
 - Single stool: 15 lbs (6.8 kg)
 - Two stools & case: 60 lbs (27 kg)
- Manufactured in the USA

NSN: 6520-01-456-5406



NEW!

ALU-40LED-CF

- Superior illumination in remote field environments
- Height-adjustable stand with gooseneck
- Adjustable high Intensity LED light with seven beam diameters
- Replaceable protective lens cover
- Amber filter for working with composites
- Super long life & durable LED
- 1,500 FC @ 24" (16,145 LUX @ 61 cm)
- Universal voltage input (100-240 VAC 50-60 Hz)
- Auxiliary DC input (12-24 VDC)
- Weight: 12 lbs (5.4 kg)
- Folded Dimensions: 9" x 45" x 5.5" (23 x 114 x 14 cm)
- 12-24 VDC Vehicle lighter plug adaptor & Battery Clips
- Includes U.S. (110V) & Euro (220V) power cords

ADC-01CS Portable Dental Chair

Professionals involved in portable dentistry continuously praise this chair! Lightweight and very sturdy. Comfortable and safe for up to 500 lbs (225 kg). High-strength steel construction. Backrest adjusts, with patient in chair, from full upright to supine position. Adjustable seat height from 11" to 22" (28 - 56 cm). Can be used as standard sit-up chair.

- High-strength DOM steel frame
- Corrosion-resistant stainless hardware
- Quality vinyl covering
- Easy, quick set-up
- Adjustable pillow headrest
- Includes hard & soft cases
- Manufactured in the USA
- Dimensions: 39" x 30" x 16.5", 11.2 cu. ft. (99 x 76 x 42 cm)
- Weight:
 - Single chair: 30 lbs (13.6 kg)
 - In mil-spec. case: 80 lbs (36 kg)

NSN: 6520-01-446-3783



ATC-03CF Portable Tray Stand

- Tripod stand with standard 13" x 9.25" stainless instrument tray on top (33 x 24.8 cm)
- Removable tray for autoclaving
- Stand is height adjustable from 29" to 57" (73.7 to 144.8 cm)
- Folds-up in vinyl carrying case. Fits in ADC-01CS mil-spec. hard chair case
- Manufactured in the USA

NSN: 6520-01-446-6263



NOTE: ADC-01CS Portable Dental Chair includes the hard mil-spec. case.

Optional items which fit into the case are as follows:

- ATC-03CF Tray stand in vinyl bag
- ALU-29CF Light in soft carrying case
- ARU-01S Tripod stand for X-ray in vinyl bag
- Combined weight: 103 lbs (46.4 kg)



Aseptico Portable Sinks

- Hands-free washing
- Can be operated electrically or with manual foot pump
- Integrated electric water heater
- Stainless steel side tray mounts to either side 6.5" x 10" (16.5 x 25.4 cm)
- Use standard military water containers for water supply and waste disposal
- Water and waste lines designed to protect against cross contamination
- Water spigot raises and lowers for ease of storage
- Legs quickly snap in and out of base and store in case lid
- Designed for easy maintenance
- Water flow approximately 2.2 liters per minute
- Dual-voltage, auto-switching; 120/220 V - 50/60 Hz.
- Manufactured in the USA



ADU-40CF Portable Field Sink

- Dimensions: 19" x 19" x 15", 3.13 cu. ft. (48 x 48 x 38 cm)
- Weight: 50 lbs (22.5 kg)

NSN: 6520-01-456-5380

ADU-40CF Optional Accessories

- **AA-04CF**
Leg brace kit NSN: 6520-01-537-4224
- **AA-04SD**
Soap dispenser NSN: 6530-01-537-7355
- **AA-04EXS**
Extended spout NSN: 6520-01-534-7852



ADU-42CF Portable Scrub Sink

- Large stainless steel basin 21.3" x 15.4" x 7.4" (54 x 39 x 19 cm)
- Stainless steel strainer basket for instruments
- Foot pump activated soap dispenser
- Paper towel holder
- Stainless steel braced legs
- Dimensions 25.25" x 22.375" x 12.625" (64 x 57 x 32 cm)
- Deck height: 32.25" (82 cm)
- Tank volume: 1.5 gal. (5.7L)
- Storage weight (incl. case): 97 lbs (44 kg)
- Weight assembled: 50 lbs (23 kg)

NSN: 6530-01-572-6775



Aseptico Products Are Proudly Engineered And Manufactured In The USA

Aseptico Military Field Dental Equipment is Available on ECAT at www.dmmonline.com



AEU-425CF Transport MDS

Aseptico's lightweight Transport MDS Electric Military Field Dental System is a result of years of analysis, design, engineering and testing. The Transport MDS provides a broad range of dental procedures to be performed in the field including General Dentistry, Oral Surgery, Endodontics and Dental Hygiene.

- 3-way air/water syringe
- Autoclavable 30,000 rpm micro-motor with internal irrigation
- Fiber optic micromotor features internal irrigation
- Accepts all E-type handpieces
- High & low volume evacuation
- Self contained water supply
- Self contained air compressor
- AE-7P On/Off Footswitch
- High impact case
- Four ratio settings: 1:5, 1:1, 4:1, 16:1
- Piezo ultrasonic scaler
- Dual voltage for world-wide operation
- UL, CUL, CE certified
- **NSN: 6520-01-508-9583**
- **Optional AE-42INV** power inverter with NATO connector allows the Transport MDS to be operated from 24 volt military vehicle's power generation
- **Optional version AEU-425CFH** includes five handpieces

Transport MDS

Transport MDS with Inverter & Mil. Spec. Case

Dimensions:	10.5" x 14.5" x 24"
Volume:	1.81 cu.ft.
Weight:	52.5 lbs (23.6 kg)
Dimensions:	25.5" x 24" x 12"
	(63.5 x 60.96 x 30.48 cm)
Volume:	4.25 cu.ft.
Weight:	80.1 lbs (36 kg)
Power Inverter:	24 VDC to 120 VAC

Power Rating:	800 W
Operating Pressure:	45 to 50 PSI
Constant Air Flow Pressure:	40 PSI, 0.4 SCFM
High Volume Vacuum:	3.2 SCFM@4.0" Hg
Low Volume Vacuum:	1.3 SCFM@1.5" Hg
Air Reserve Capacity:	61.0 cu. in (1000 cc)
Noise Level, "A" Weighting:	70 dB@3'4" (1m)



Model AEU-14CF

In response to the Military's demand for an ultra lightweight emergency dental unit, Aseptico has developed the AEU-14CF Expedition Field Dental System.

Weighing under 25 pounds, the Expedition is light and compact enough to be carried in a backpack and deployed into the field of operation along with the medical support team.

This unique electric dental system can be powered via AC power source, 12 or 24 V vehicle battery, or from its two rechargeable battery packs. For maximum flexibility in the field, the battery packs may be recharged using the included solar panel, vehicle battery, or AC power source. Each battery pack provides power for approximately 20 patient procedures.

- AHP-72 1:5 Increaser Handpiece
- AHP-101 1:1 Straight Handpiece
- Handpiece irrigation system
- Autoclavable 30,000 rpm micro-motor with internal irrigation
- Fiber optic micromotor features internal irrigation
- 3-way air/water syringe
- Integrated compressor for irrigation & syringe
- Operates from rechargeable battery packs, AC outlets, 12 or 24 V vehicle battery
- Battery packs recharged via solar panel (*included*), vehicle battery, or AC power source
- Dual voltage for world-wide operation
- Weighs under 25 lbs (11.25 kg)
- Dimensions: 16.2" x 12.7" x 6.5" (41 x 32 x 16.5 cm)
- Noise Level: 64 dB
- Water Reservoir: 16 fl. oz.
- Air Reservoir: 16 fl. oz.
- UL, CUL, CE certified

AMC-20CF Mobile Dental Cart

- Fully self-contained dental cart with built-in air compressor/vacuum pump, plus water & waste systems
- 2 fiber optic highspeed handpiece connections
- Integrated AEU-5000 Electric dental motor (see page 3)
- Powerful suction with overflow protection, electric waste purge & amalgam separator
- Piezo ultrasonic scaler
- Curing light
- Mil-Spec case with ramp
- Large countertop with height adjustable tray & pivoting control arms
- Dual voltage for world-wide operation
- UL 60601-1 certification for medical device safety compliance (*pending*)
- Dimensions:
L 30" x W 23.5" x H 36.5"
(76.2 x 59.7 x 92.7 cm)
- Case Dimensions:
L 45" x 32" x 30.4" (114.3 x 81.3 x 77.1 cm)
- Weight: 194 lbs (88 kg)
- Case Weight: 210 lbs (95.3 kg)
- Power Rating:
115 V - 1,620W @ 60 Hz or 1,944W @ 50Hz
230 V - 1,725W @ 60 Hz or 1,437.5W @ 50Hz

**ARU-07LC NOMAD™ PRO with Large Case****ARU-07 NOMAD™ PRO with Standard Case****ARU-06CEL NOMAD™ CLASSIC with Large Case****ARU-06CE NOMAD™ CLASSIC with Standard Case**

These hand-held X-ray systems are powered via rechargeable battery packs, similar to those used in cordless drills.

The X-rays are lightweight and have internal and external shielding which allows the operator to safely take the exposure while standing with the patient.

See page 16 for features and accessories. Please visit us at www.aseptico.com for complete specifications and manuals.

ARU-06CE NOMAD™ CLASSIC
NSN: 6525-01-551-3196



NOMAD™ PRO



NOMAD™ CLASSIC



STANDARD CASE

EXCLUSIVELY FROM ASEPTICO!

Optional large case with space for NOMAD™ PRO or NOMAD™ CLASSIC, laptop, sensors & all accessories



Piezoelectric Scaler



Replacement Parts

AA-12HP Replacement Handpiece	\$250
AE-7PM On/Off Footswitch.	\$125

ASC-12 Piezoelectric Scaler \$995

ASC-12 with 3 Tips. \$1,200

- Wide variety of uses, including: Hygiene, Periodontal & Endodontic scaling applications
- 16 Power levels (5 power boost settings)
- Programmable preset buttons
- Compatible with a wide assortment of scaler tips (M3 x 0.6 thread pattern)
- Auto-sensing global voltage compatibility
- AE-7PM 360° touch, weighted footswitch
- Auto-detects frequency range: 27 - 34 kHz
- Manufactured in the USA

	American Eagle Infinitip™	Equivalent Satelec® Tip	Suggested Power Level										
				Light - Medium Calculus	High Density Calculus	Furcation	Shallow Pockets	Supra Gingival	Sub Gingival	Molars & Pre-Molars	All Tooth Surfaces	Stain Removal	
Scaling	ASC-12-11A	1	Med	✓				✓					
	ASC-12-12B	2	Med		✓			✓					
	ASC-12-13C	3	Med								✓	✓	
	ASC-12-14P	10P	Med			✓							
	ASC-12-14U	PFU	Med					✓					
	ASC-12-14R	PFR	Med					✓					
	ASC-12-14L	PFL	Med					✓					
	ASC-12-14X	10X	Med	✓				✓					
	ASC-12-14Z	10Z	Med	✓				✓					
Periodontics	ASC-12-14D	H1	Low	✓									
	ASC-12-14BLD	H2L	Low		✓				✓				
	ASC-12-14BRD	H2R	Low		✓				✓				
	ASC-12-14DL	H4L	Low		✓				✓	✓			
	ASC-12-14DR	H4R	Low		✓				✓	✓			
	ASC-12-14C	H3	Low						✓				

American Eagle Infinitip™ XP Technology

- No loss of efficiency - the last stroke feels like the first
- Designed to be the longest-lasting tip on the market
- Improves performance and extends the life of control units
- Offers greater patient and user comfort
- Satelec® Tips for reference only - not available through Aseptico
- AEI Infinitips™ made by American Eagle Instruments®
- Available to US customers only
- Satelec® and Infinitip™ are trademarks of their respective owners



ASC-12-14P
\$95



ASC-12-14U
\$95



ASC-12-12B
\$95



ASC-12-14R
\$95



ASC-12-13C
\$95



ASC-12-14L
\$95

RECOMMENDED STARTER TIPS



ASC-12-14Z
\$95



ASC-12-11A
\$95



ASC-12-14C
\$95




ASC-12-14X
\$95



ASC-12-14DL
\$95




ASC-12-14DR
\$95

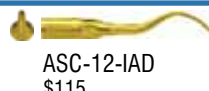

ASC-12-14D
\$115



ASC-12-14BRD
\$115

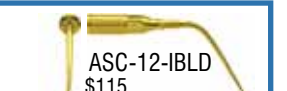

ASC-12-14BLD
\$115

DIAMOND COATED TIPS